

Xialla

What is Xialla?


Xialla is a loop tension band system that helps to maintain an erection when there is a venous leak. Venous leak is a condition where the penis cannot keep blood inside of it to create or sustain an erection. The Xialla consists of four parts: The ring, loop, band, and tab lock. An adhesive pad is also required and comes with the product. All parts work to secure Xialla, so that it stays in place during sexual play.

You can find the product and sizing information on the website:


ca.xialla.com/products/xialla

Available sizes include: Regular, Slender, and Thicker.


Discount code: SG55683 (10% off)

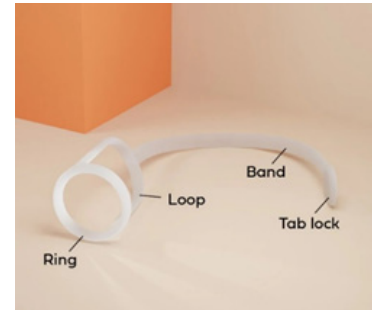


How to Wear Xialla

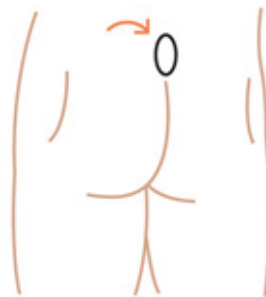


1. Clean and dry your lower back. Remove protective backing from the adhesive pad and stick the pad above your buttocks.

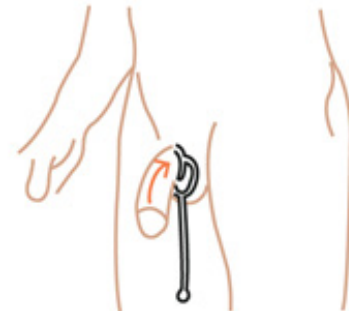

2. Place the ring around the base of the penis. Let the loop and band hang in front of the scrotum. The ring can be placed on an erect or flaccid penis.

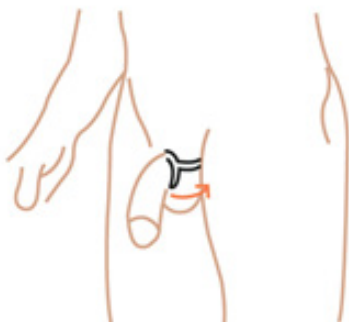

3. Stretch the loop up and around the scrotum so that it is inside the loop. Adjust the loop to fit comfortably above the scrotum.

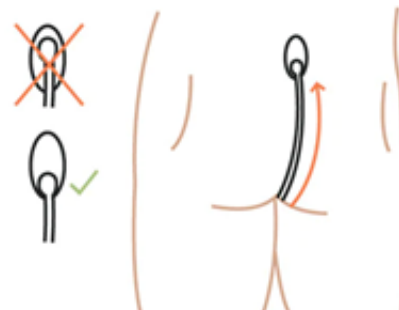

4. Pull the band back between your legs and attach the tab lock to the bottom of the adhesive pad to lock Xialla in place.

Video on how to apply the Xialla: youtu.be/bwdntlpMgYw

Can Xialla be used with a vacuum erection device (VED)?

Yes, Xialla can be used with VEDs. Please see instructions below.

For the VED where the tension band is applied over the cone:

1. Apply the adhesive pad to your lower back.
2. Place Xialla ring on the pump using the cone.
3. Pump the VED as you normally would.
4. Slip the Xialla ring off the pump onto the penis and remove the pump.
5. Finish putting on the Xialla by placing the loop around the scrotum, pulling the band back between your legs, and attaching the tab lock to the adhesive pad.

For the VED with a tension ring that forms the vacuum seal:

1. Apply the adhesive pad to your lower back.
2. Pump as you normally would and then remove the VED.
3. Place the Xialla ring between the VED's ring and the base of the penis.
4. Finish putting on the Xialla by placing the loop around the scrotum, pulling the band back between your legs, and attaching the tab lock to the adhesive pad.
5. Remove the ring that came with the VED.

**Please note, all tension bands must be removed after 30 minutes to allow blood to circulate.*

Images used with permission from Xialla

For more copies, go online at vch.eduhealth.ca or email pchem@vch.ca and quote Catalogue No. **GD.352.X53**

© Vancouver Coastal Health, June 2023

The information in this document is intended solely for the person to whom it was given by the health care team.

vch.ca



Making better
decisions together
with patients
and families